

# Testing the classifier

## iTongueO project

[Preparing the environment](#)

[Electrodes](#)


[Online classification](#)

### Preparing the environment

1. Download the project from here: [HyperDecahedron/tongue\\_emg\\_xr\\_project: Project to identify tongue's position as an XR input.](#)
2. Prepare your environment to execute a python code in tongue\_project/scripts.

### Electrodes

1. Take 7 electrodes. If the brand is Ambu, they don't require additional gel. If the brand is 3M, add a bit of conductive gel to the center.
2. Place the electrodes. The image below should be interpreted as a mirror (blue dots are your left, gray dots are your right).

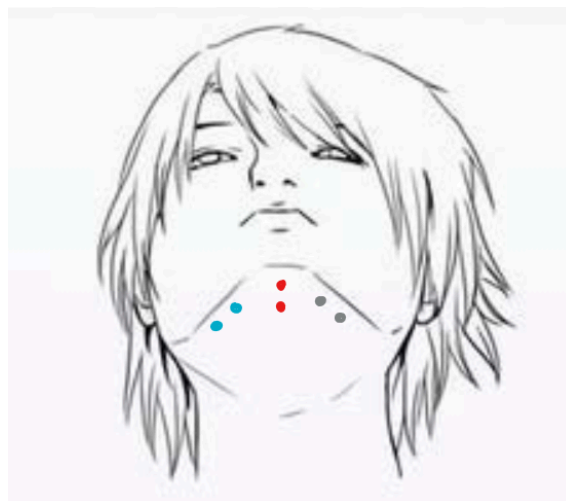


- a. Place 1 electrode (bias) in your mastoid. The wire should be black.
- b. Place 2 electrodes for **Channel 1 (red in the picture)** in the **center** of the submental triangle. If the board was not modified, the wires should be **red** and **brown** (N1P).
- c. Place 2 electrodes for **Channel 2 (gray in the picture)** in your **right side** of the submental triangle. If the board was not modified, the wires should be white and **gray** (N2P).

- d. Place 2 electrodes for **Channel 3 (blue in the picture)** in your **left side** of the submental triangle. If the board was not modified, the wires should be **blue** and **purple** (N3P).

## Online classification

1. Turn on the board, connect the dongle and open OpenBCI GUI.
2. **Set off all channels except for channels 1, 2 and 3.** To do this, click on the icon of each number (I cannot open the GUI without the board, but these little numbers should be gray from 4 onwards).

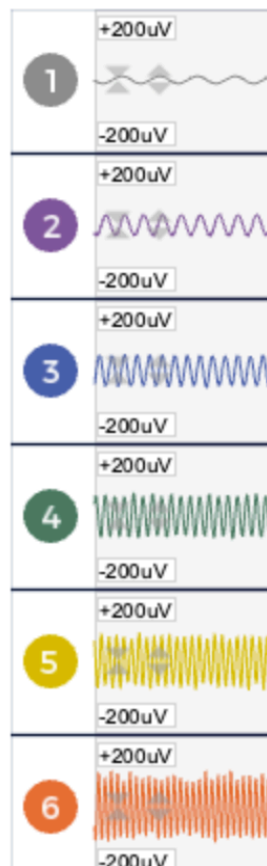

3. Open Hardware Settings and **set SRB2 Off for channels 1, 2 and 3.** Click 'Send' to save the settings
4. Open the widget for Networking and select **UDP** and **TimeSeriesRaw**.
5. Start the session.
6. Start the UDP networking.
7. Launch the python script "**real\_time\_6\_classes.py**" in **tongue\_project/scripts**.
8. Test the classes by pressing your teeth with the tip of the tongue, according to the following table.

Class	Meaning	Tooth
l	left	Upper Left Third Molar
lf	left front	Upper Left First Premolar
f	front	Upper Between the Central Incisors
rf	right front	Upper Right First Premolar
r	right	Upper Right Third Molar
s	swallow	Swallow or fake the movement by pressing your submental triangle firmly upwards
n	none	Tongue is resting

## Adult dental chart

Eruption (year)

